


BREAKAWAY BED CABLE COMPARISON: CURBELL VS. MAGNETIC

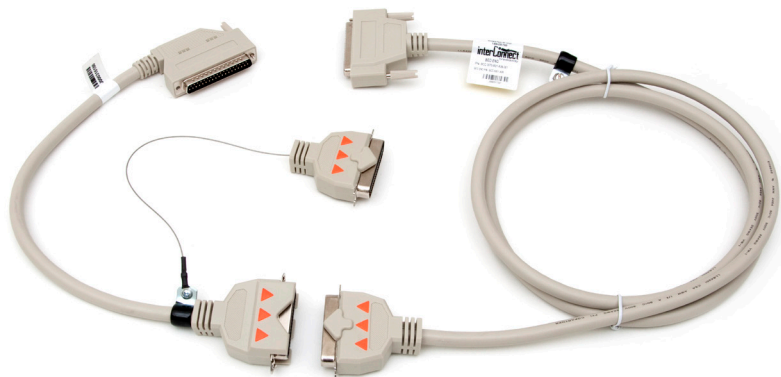

The safety of your patients is not something you want to gamble with. Curbell has used our 50 years of nurse call system expertise to design a breakaway cable specifically with safety in mind. This cable includes several features that help ensure a proper and safe connection between your patients and staff, all at a fraction of the cost of the competition.

	Curbell's Breakaway Bed Cable	Magnetic Breakaway Bed Cable	
Easy Installation	●	●	<p>Curbell's breakaway cables have one end permanently mounted to the wall and the other end permanently mounted to the bed, and the two sections connect via a universal pinless connection.</p> <p>Magnetic cables have one end permanently mounted to the bed, while the other end magnetically mounts to the wall (or vice versa) using an adaptor.</p>
Proper Disconnect	●		<p>Curbell's breakaway cables separate at the middle, with permanent mounts securely in place at the wall and the bed to eliminate unwanted disconnects.</p> <p>Magnetic cables are only permanently mounted to one end and can easily separate if someone bumps into the end that's not fastened.</p>
Price	●		<p>Magnetic cables are typically three times more expensive than Curbell's breakaway cables.</p>
Voltage/Patient Safety	●		<p>Curbell's breakaway cables have molded ends to ensure that the cable is not connected backwards, and the pinless connection won't attract metal objects that could cause a short.</p> <p>Magnetic cables can attach to metal surfaces or attract loose objects (staples, needles, paperclips, etc.), which can cause a short circuit in the nurse call system or bed. These objects can also be a source of contamination.</p>
Nurse Call System Compatibility	●		<p>Curbell's breakaway cables are available with connections for the majority of nurse call systems.</p> <p>Magnetic Cables are only available with a 37 pin or 1/4" jack connections.</p>
Customizability	●		<p>Curbell's breakaway cables are available in lengths from 1' to 15' for either the wall or bed section to customize them for your individual patient room layouts.</p> <p>Magnetic cables are only available in 1' and 6.5' lengths.</p>
Proven Reliability	●		<p>Curbell's breakaway cables have been installed in over 1000 facilities in the last 15 years.</p>