


CURBELL MEDICAL MATTRESS SOLUTIONS

SEALED SEAM MATTRESSES

Our sealed seam bed mattresses feature a fully sealed cover suitable for behavioral health applications. They feature three layers of medical-grade foam to provide a high level of comfort at a lower price than the OEM. Firm foam sides give extra support when a patient is sitting on the edge of the bed, and the standard fire barrier meets CAL129/FED1633 fire rating.

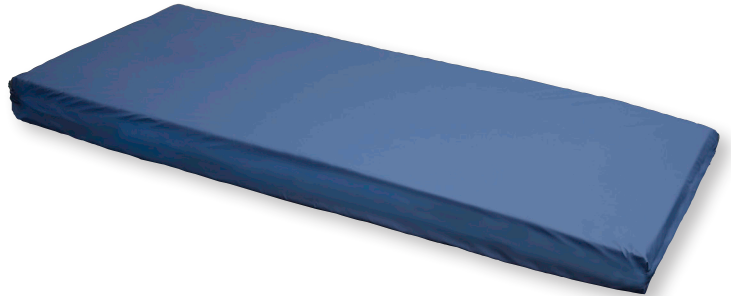

RF WELDED SEAMS

Rather than a sewn seam, our seams are sealed using radio frequency welding to prevent fluid ingress. (Only available on polycarbonate material.)



PERFECT FIT

All our mattresses are accurately made to be a perfect fit on your beds.



POLYCARBONATE FABRIC

Our mattresses feature a durable polyester-coated knit with a polycarbonate additive. This material is resistant to bleach and peroxide.



FIRE SAFETY

All our foam meets Fed 1633 fire safety standards.



BI-DIRECTIONAL STRETCH

Our breathable, two-way stretch fabric gives added comfort and durability, and also prevents “hammocking,” allowing for easier patient transfer.



FIVE YEAR WARRANTY

We stand behind all our bed mattresses with a five year warranty.



THOROUGH TESTING


All of our mattresses are tested by Support Surface Standards Initiative (S3I). Visit <https://npiap.com/page/S3I> for more information on this organization.



MADE IN THE USA

All of our mattresses are manufactured in the USA.

STANDARD FOAM CONSTRUCTION

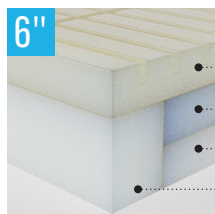

- 2" Conventional Foam (sculpted cubes)
- 2" High Resiliency Foam
- 2" Conventional Foam
- 4" Conventional Foam Side Rail

				FULL MATTRESS	
Width	Length	Thickness	Patient Weight Limit	Part Number	
35"	80"	6"		28748	
36"	80"	6"		28743	

©2024 Curbell Medical Products, Inc MAP1875C