





EMPOWERS PATIENTS AND CLINICIANS

Patients are desiring greater social connectivity and immediate access to information. As a result, some hospitals are providing commercial tablets to patients to enhance the hospital stay. These consumer tablets introduce challenges with connectivity, charging, device tracking and management, as well as data security.

As the leader in patient communication devices for over 50 years, Curbell Medical proudly introduces the new Rēgo™ Patient Interaction System — a scalable platform that addresses the challenges associated with tablets, integrating all patient controls into a single, medical-grade device.

The Rēgo™ system is designed to meet patient and hospital requirements with a focus on enhanced communication, information access, security, durability, cleanability, and regulatory compliance.

RēGO IS PURPOSE-BUILT FOR THE PATIENT ROOM



Rēgo is an Android™ device. Curbell is licensed for Google Mobile Services and certified as a Google Original Device Manufacturer.



Rēgo is a medical device per the FDA and meets all compliance requirements including UL 1069 and part 15 of the FCC rules.



Rēgo's display is a 7" diagonal touchscreen with a resolution of 1024 x 600 and is protected by Corning® Gorilla® Glass.



Designed to survive 100 drops from four feet onto a concrete floor and compatible with most common hospital cleaners and disinfectants.



Network security and patient privacy is paramount. You decide what apps to install, and wiping data between patients is quick and easy.



Bluetooth functionality provides the ability to connect Rēgo with up to 10 additional devices in the room.



Rēgo can connect to your hospital's network via WiFi or hardwired ethernet.



Rēgo includes a front facing camera for video chat. A 5 MP rear-facing camera with flash is optional.





NURSE CALL IS THE CENTER OF THE PATIENT EXPERIENCE

Curbell has been making products to integrate the systems in the patient's room since 1962. By combining our unparalleled institutional knowledge and decades of experience, we create innovative solutions that make hospitals more efficient, and improve patient comfort and safety. We know that empowering patients starts with giving them the peace of mind knowing that they can quickly and easily call the nurse via your nurse call system, and with Rêgo, calling the nurse is just the beginning of what patients can do.

RĒGO IS CHANGING THE PATIENT EXPERIENCE

Practical and purpose-built for hospital patients and staff, Rēgo is easy to see and use for patients of all ages. Rēgo empowers patients to communicate their needs, select entertainment options, access information, and even control in-room elements like lighting and shades.



The Rēgo App Marketplace offers apps curated specifically to meet the needs of hospitals and patients — with privacy and security in mind. Combined with the Rēgo hardware, hospitals are now dramatically enhancing the patient experience and increasing clinical efficiency.



INTERNET

With Rēgo's built in web browser, patients can access the internet to stay connected via social media, news, weather, etc.



GAMES

Provide another valuable entertainment alternative for patients — ideal for passing the time and reducing anxiety.



PATIENT ENGAGEMENT SYSTEMS

Rēgo can integrate with your existing patient engagement system, giving patients touchscreen access to the system's features in addition to the TV. Control it all from a single device!





ELECTRONIC MEDICAL RECORDS

Leverage Rêgo to integrate with your existing EMR system to provide valuable information to patients and increase clinical efficiency for your staff.



E-BOOKS

Electronic books can provide an alternative entertainment option for patients.



EDUCATION

Give your patients access to documents and videos that help them better understand their condition and treatment.



SURVEYS

Gather feedback from patients before discharge, while you can react and ensure a positive experience.



THE TELEVISION EXPERIENCE

Patients control the TV with the dedicated channel, volume, mute, and power buttons, while pop-up screens provide control of additional television features. Audio comes through Rêgo's stereo speakers, but is silenced during nurse call communications.






VIDEO CHAT

Patients can connect via video chat with friends and family. Your hospital can even use the system to allow patients to communicate with doctors or a translation service.



ROOM COMFORT CONTROL

Give your patients enhanced control of room lighting, shades, temperature, and more from the dedicated room control screen, customized to the room level.

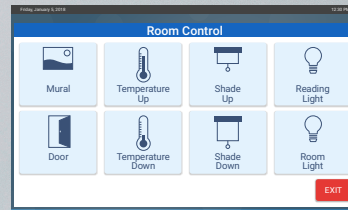

TRANSLATION

Leverage your existing translation service to be delivered via Rêgo in every patient room — always available when needed without additional equipment.



All buttons are back-lit, making them easy to find in the dark. Additional LEDs give visual assurance of various functions.




The device is powered and tethered at all times. When not in use, it can be placed into the holder, which can mount to most bed side rails.



The same holder can be used as a stand to allow use on an over-bed table.

THE RÊGO SYSTEM INSTALLS IN MINUTES

Installation can be accomplished without room downtime and does not require access behind the wall.

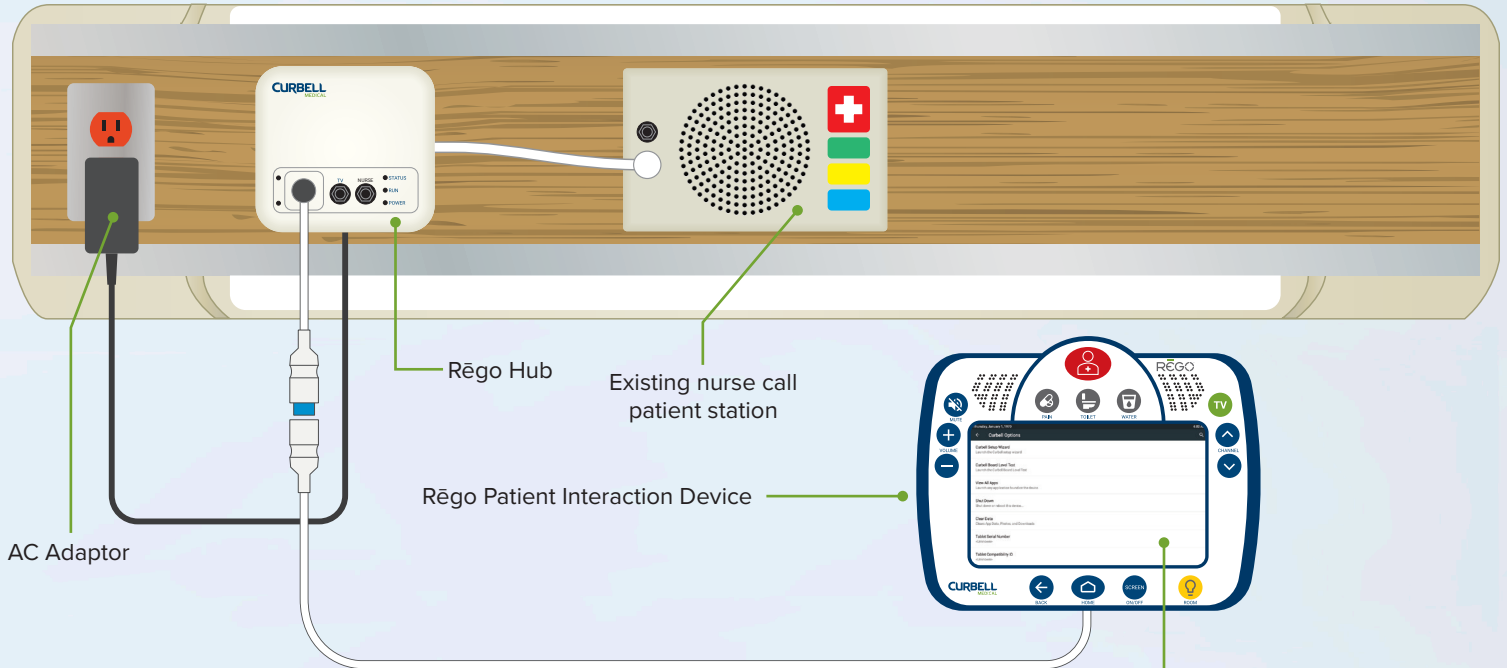


SETUP AND SECURITY

Configuring the device for your network is simple, using the Setup Wizard. When a patient is discharged, all user data is easily deleted, ensuring that anything entered into the device during the patient's stay is deleted upon discharge.

THE RÊGO HUB

The Hub is a key component to the Rêgo system that performs several functions:

- 1 Provides constant power to the Patient Interaction Device (PID) so it never needs charging
- 2 Connects to your existing nurse call system
- 3 Provides connections for additional hardwired room controls (lights, shades, etc.) to supplement what the nurse call station offers
- 4 Provides an ethernet connection option in addition to the PID's built-in WiFi

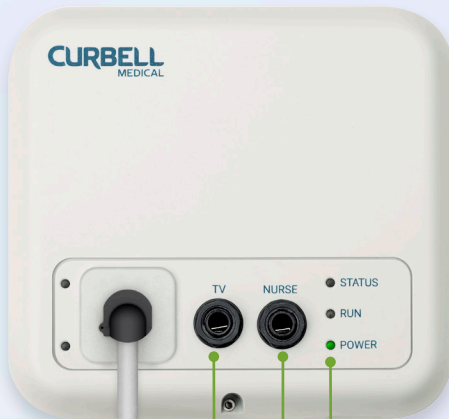

LEDs provide visual verification of system status

1/4" jacks accommodate specialty call devices

Rêgo's breakaway connector reduces plug damage issues

THE RĒGO ECOSYSTEM ENSURES PEACE OF MIND

RĒgo is purpose-built for hospital patient rooms and designed to ensure security and patient privacy. The RĒgo App Marketplace, combined with the on-board Mobile Device Management agent, makes it easy for your IT administrator to manage devices remotely from a single computer.




RĒGO APP MARKETPLACE

Connect to Curbell's RĒgo App Marketplace to choose from a wide selection of apps designed for use by patients during their hospital stay. An SDK is available for app developers to leverage the platform.



SYSTEM ADMINISTRATION

The hospital downloads apps and software updates from the marketplace to the Patient Interaction Devices via a Mobile Device Management system to ensure safety and security.

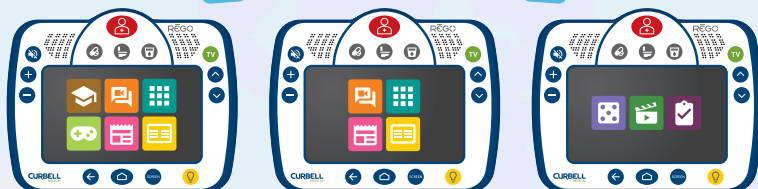

REMOTE VISIBILITY AND CONTROL

The hospital maintains visibility and control of units to the patient room level, with the ability to offer specific apps to individual patients.

Call today to schedule a personalized webinar to learn more.

CURBELL
MEDICAL

THE (R)EVOLUTION OF PATIENT CONTROL DEVICES



In hospital rooms around the world, patients are dealing with demands on their health, and staff members are dealing with demands on their time. Curbell Medical is at the center of it all, with product lines that have set the standard for hospital room integration. For more than 50 years, we have made sure that patients can expect the best possible experience.

www.hellorego.com • 1-800-235-7500